

icanCLAVE



ICANCLAVE E series  
Class N autoclave 18L, 23L



### Sterilization cycles

Two sterilization programs, 121°C/1.1Bar and 134°C/2.1Bar.



### Time

Adjustable sterilization and drying time.



### Pressure safety lock

It only allows the door to open when the pressure is 0.



### Drying (option)

Heat drying function.



### Automatic pressure release

Automatically discharge of water and steam after cycle over.



### High precision control

High precision temperature control within 0.1 °C.



### Protection

Temperature and pressure overload protection.



### Safety door lock

It only allows the cycle to start when the door is properly locked.



### Drying + (option)

Drying plus function for reducing drying time and better drying result.



Icanclave E is class N steam sterilizer series which includes the chamber size of 18L, 23L, and 50L, they provide an economic sterilization solution for dental clinics, laboratories, cosmetic hospitals, tattoo shops, veterinarians, etc.

E - X X X

18: 18Liter; 23: 23Liter;

3: no drying 4: with drying\* 5: with drying and drying plus\*

\*drying: drying by heater. For better drying results, open the door and leave a small gap once the autoclave finishes the cycle.

\*drying plus: it equips with an additional air compressor that is activated automatically once the cycle is finished improving drying results.

Fox Example : E418 means E series 18L autoclave with drying function

## How to order

# E\*18, E\*23

## Technical Parameter

Order Code	E318/E323	E418/E423	E518/E523
Model Number	STE-18-E(18L)/ STE-23-E(23L)		
Chamber Size	Φ250x350mm(18L)/Φ250x450mm(23L)		
Electrical Supply	220~240 V/50Hz, 110~120V/60Hz		
Power	1500VA		
Sterilization Temp.	121°C / 134°C		
Sterilization Time 121°C	Default 30min (30~60min adjustable)		
Sterilization Time 134°C	Default 10min (10~20min adjustable)		
Pressure	1.1Bar/121°C, 2.1Bar /134°C		
Automatic Discharge	Yes	Yes	Yes
Dry Time	No	Default 60min (1~60min adjustable)	
Dry+ Time	No	No	Default 60min 1~60min adjustable
Net Weight	32 kg(18L)/36kg(23L)	33 kg(18L)/37kg(23L)	33.5 kg(18L)/37.5kg(23L)



Heater



Safety Lock



Operating Membrane




Discharge Port

Power Switch


Display Window

# Installation & Overall Dimension

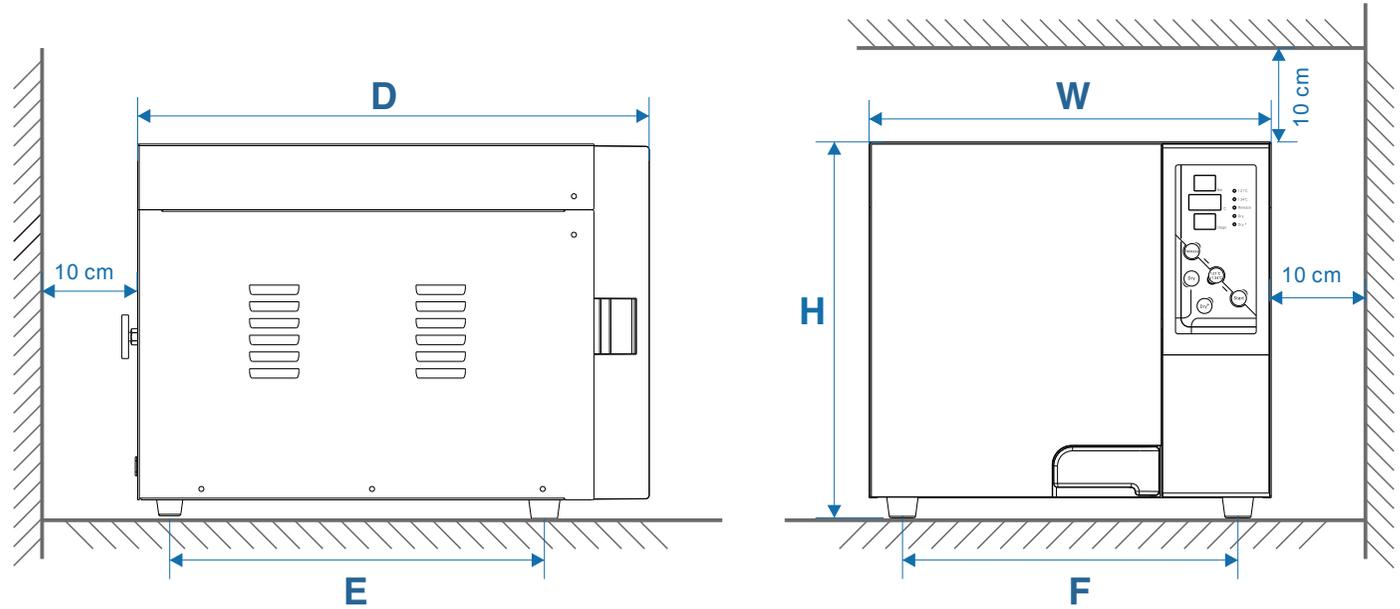


Item	W	H	D	E	F
E*18	483mm	450mm	515mm	350mm	403mm
E*23	483mm	450mm	615mm	450mm	403mm

# Rack & Tray Dimension

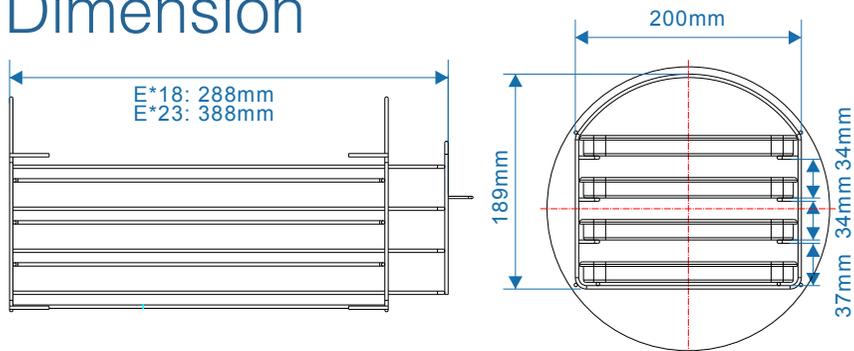

# Tray Dimension

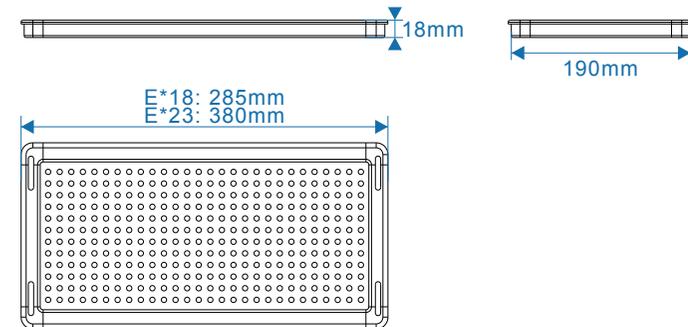

\* Number of Trays: 2 Pcs Standard, 4 pcs Option.

